PRODUKTNAME
PRODUCT NAME

System SIN

BEZEICHNUNG
DESCRIPTION

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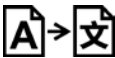
Wichtige Hinweise / Important information

Intended use

System SIN is a medical device for the additive manufacturing of dental frameworks using the laser melting process.

Only for professional user (Dental Technician, Dentist).

The intended patient group provides for persons with partially or non-dentate jaws.

Indication

Crowns and bridges for metal-ceramics, frameworks for removable dentures, superstructures and bars.

Contraindication

In case of known allergic reactions to any of the ingredients.

CAD

The design should be done with appropriate CAD software. Please consider an anatomically reduced framework design for the veneering with ceramic. The wall thickness should not be less than 0.3 mm. Choose a sufficient connector dimension (6 - 9 mm²). Sharp edges and undercuts should be avoided.

Thermal treatment

After the laser melting process, the building boards have to pass heat treatment to minimize stresses due to the laser sintering process. For this purpose, a suitable furnace with inert gas (argon) or vacuum function should be used. Stress relieving without inert gas atmosphere can optionally be performed. Please consider that a treatment without inert gas atmosphere can lead to an increased oxide formation.

Description	Temperature [°C]	Time [min]
Stress-relief-heat-treatment with inert gas	850	30
Cooling 1 with inert gas	Cooling to 600, then open door	
Cooling 2, w/o inert gas	300 – room temperature	

Remove Parts from the build plate

After heat treatment and cooling the plate, the restorations can be removed by band saw, rotary instruments or piers.

Soldering / Laser welding

System SIN can be soldered with all suitable solder. **System SIN** parts should not be soldered with gold or palladium solder. **System SIN** is also ideally suitable for laser-welding.

Preparation before ceramic veneering

The frameworks can be elaborated with standard carbide cutters, look for smooth transitions and avoid overlapping material. Please use the same cutter for one alloy to avoid contamination. The minimum thickness of the prepared coping should not be less than 0.3 mm. It's recommended to sandblast the frames with minimum 110 µm of Aluminium oxide with 3-4 bar and clean with steam cleaner. Oxide firing is not mandatory but can be done as an option for 5 minutes at 980 °C with vacuum (cleaning firing). The frame needs to be sandblasted with aluminium oxide 110µm and 3-4 bar to remove the present oxide layer thoroughly. In the end the cleaning by steam cleaner is mandatory. If you use a ceramic bonder, please consider the instruction for use of the manufacturer.

Handling conditions / Safety

Metal powder or dust may cause irritation by inhalation and in contact with skin. During the handling with the System SIN as well as while grinding and sandblasting of the produced units it's recommended to consider an adequate extraction system, goggles, gloves, protective clothing and a respirator with fine particle filter (type FFP3 – DIN EN 149). After working with the powder or restorations clean hands thoroughly.

Residual Risks and side effects

If the instructions are observed during the production processes, incompatibilities with non-precious dental alloys are extremely rare. In case of a proven allergy against an ingredient of this alloy, the alloy must not be used for safety reasons. In exceptional cases, electrochemically induced, local irritations have been reported. When different alloy groups are used, galvanic effects might occur. Please inform your dentist about residual risks and side effects. Any serious incident that involves the product must be reported to the manufacturer and the competent authority in the accorded country. The SSCP is available at <https://ec.europa.eu/tools/eudamed>.

Disinfection of the dental prosthesis before insertion

Workpieces from the dental laboratory must be subjected to immersion or spray disinfection before insertion into the patient's oral cavity and then rinsed under running water.

Disposal Instructions

Please dispose of metal residues and dust in an environmentally friendly manner. Do not allow waste to enter groundwater, water or sewage systems. Contact waste exchanges for recycling. Outer packaging can be disposed of in paper waste.

Storage conditions

Keep dry in a sealed container.

Our information and recommendation are based on the state of the art in science and technology and has to be considered correct to the best of our knowledge and experience on this day. The above version shall replace any previous versions.

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Adentatec GmbH